 FDA_Submission.md

FDA Submission

Your Name: Francois Masson

Name of your Device: X-Ray Pneumonia Detector III

Algorithm Description

1. General Information

Intended Use Statement:

For assisting a radiologist in detection of pneumonia in x-ray images.

Indications for Use:

Screening of x-ray images.

Patient population:

- Both men and women
- Age: 2 to 90

X-Ray image properties:

- Body part: Chest
- Position: AP (Anterior/Posterior) or PA (Posterior/Anterior)
- Modality: DX (Digital Radiography)


Device Limitations:

The model is recommended for use without following comorbid thoracic pathologies:


- Consolidation
- Edema
- Effusion
- Hernia

This can be explained by the fact that X-Ray images containing four conditions mentioned above have similar but significantly different distribution that that of Pneumonia X-Rays:

Pneumonia:

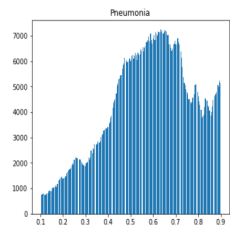

Not recommended:

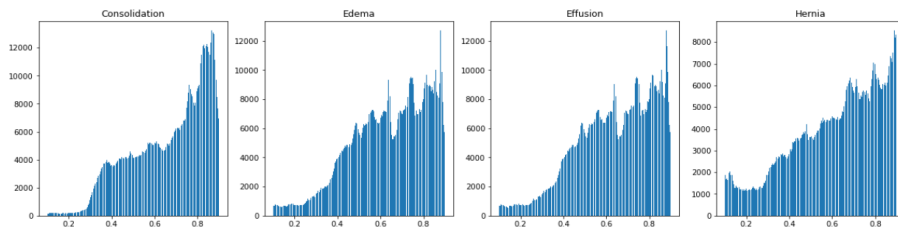

Hardware Requirements:

In terms of computer hardware requirements, the algorithm's performance was measured on Intel(R) Xeon(R) @ 2.30GHz CPU:

- Average image pre-processing time: 29 ms, max: 111 ms
- Average CNN inference time: 615 milliseconds, max: 734 ms

So, total inference time does not exceed 850 milliseconds on Intel Xeon CPU. Similar (or higher performance) CPU is recommended.

Clinical Impact of Performance:


In terms of predictive value:

- If the model predicts negative, it is correct with 90.5% probability
- If the model predicts positive, it is correct with 26.8% probability

Therefore, the algorithm is recommended for assisting a radiologist screening those images that most probably do not contain Pneumonia, to prioritize his/her time and attention to those that potentially do. This should lead to those that require medical attention getting it much faster.

It is also worth mentioning that when algorithm predicts negative is can still be wrong with 9.5% probability. So, those cases predicted negative should still be reviewed by the radiologist.

2. Algorithm Design and Function

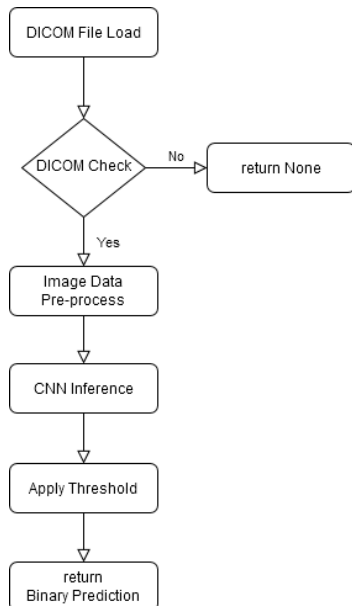

DICOM Checking Steps

The algorithm performs the following checks on the DICOM image:

- Check Patient Age is between 2 and 90 (inclusive)
- Check Examined Body Part is 'CHEST'
- Check Patient Position is either 'PA' (Posterior/Anterior) or 'AP' (Anterior/Posterior)

- Check Modality is 'DX' (Digital Radiography)

Preprocessing Steps

The algorithm performs the following preprocessing steps on an image data:


- Converts RGB to Grayscale (if needed)
- Re-sizes the image to 244 x 244 (as required by the CNN)
- Normalizes the intensity to be between 0 and 1 (from original range of 0 to 255)


CNN Architecture

The algorithm uses pre-trained VGG16 Neural Network (except the last block of Convolution + Pooling layers that was re-trained), with additional 4 blocks of 'Fully Connected + Dropout' layers.

The network output is a single probability value for binary classification.

Below is the CNN architecture graph:

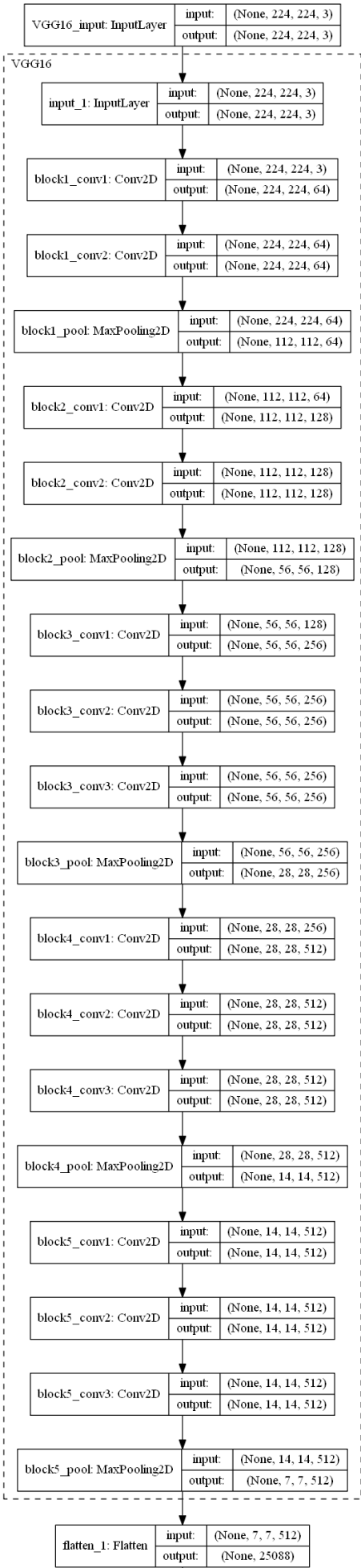

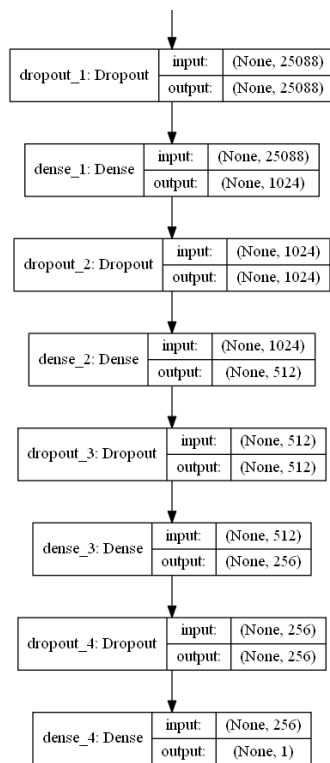


3. Algorithm Training

Parameters:

- Types of augmentation used during training:
 - horizontal flip
 - height shift: 0.1
 - width shift: 0.1
 - rotation angle range: 0 to 20 degrees
 - shear: 0.1
 - zoom: 0.1
- Batch size: 32
- Optimizer learning rate: 1e-5
- Layers of pre-existing architecture that were frozen
 - All except last convolution + pooling block
- Layers of pre-existing architecture that were fine-tuned
 - The last 2 layers of VGG16 network: block5_conv3 + block5_pool
- Layers added to pre-existing architecture
 - flatten_1 (Flatten)
 - dropout_1 (Dropout)
 - dense_1 (Dense, 1024)
 - dropout_2 (Dropout, 0.2)
 - dense_2 (Dense, 512)
 - dropout_3 (Dropout, 0.2)
 - dense_3 (Dense, 256)
 - dropout_4 (Dropout, 0.2)
 - dense_4 (Dense) 1


Algorithm training performance visualization

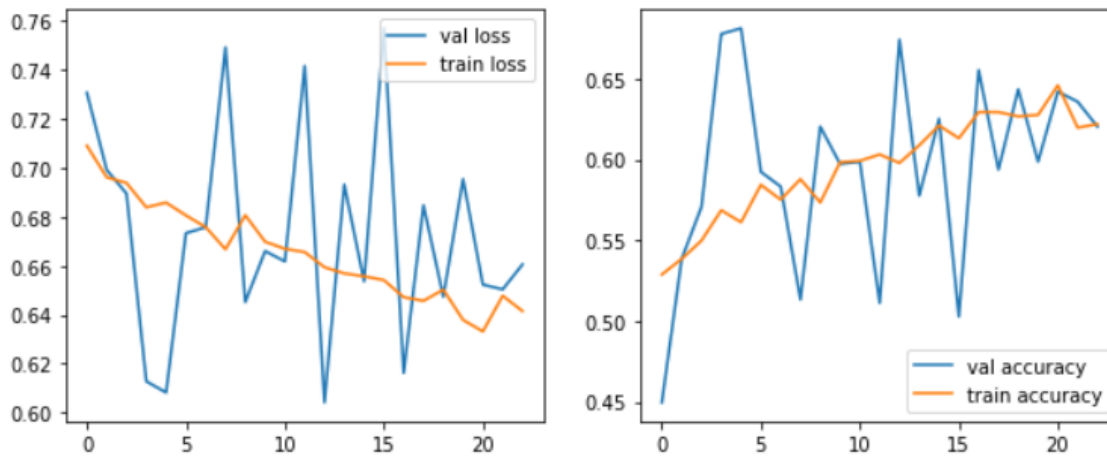


The early stopping with patience of 10 epochs was used to stop the learning process.

The behaviour of the validation loss during training may indicate a possibility that a lower learning rate could produce better results. This is left for future research.

Model performance metrics depending on threshold

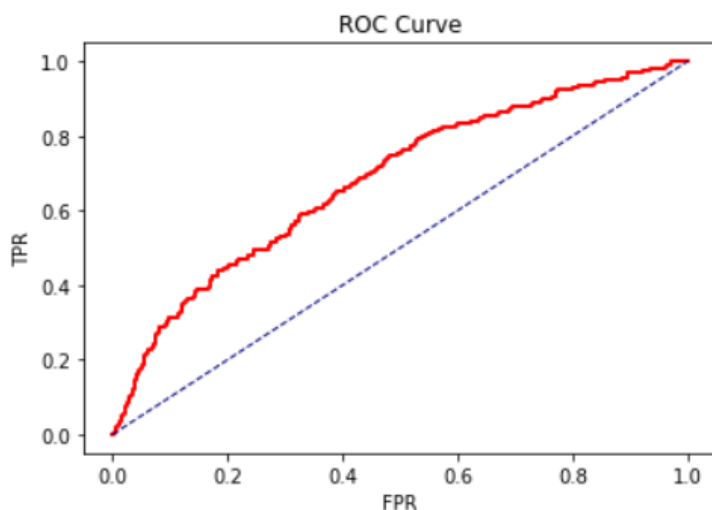

ROC Curve is not particularly impressive, but does show that the model has indeed learned something from the data.

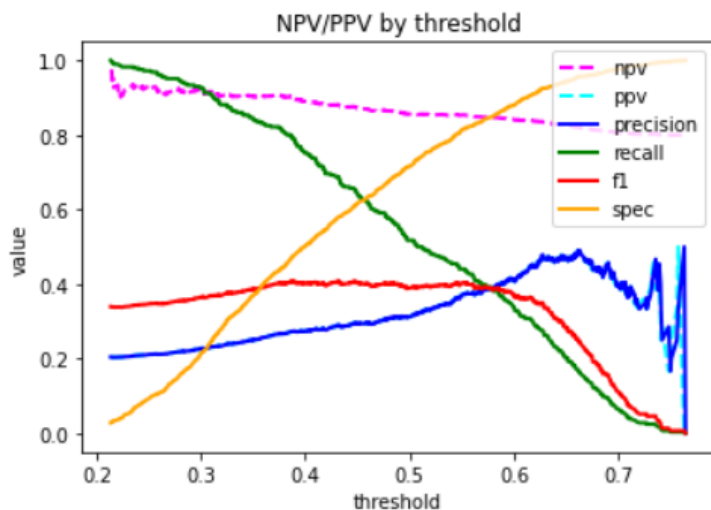

As we can see, the model has low precision, but higher recall, and maintains high negative predictive value.

Final Threshold and Explanation:

The maximum F1 score for the model is 0.408 and it is achieved with threshold value of 0.431. Below is the comparison of F1 score with those given in [CheXNet: Radiologist-Level Pneumonia Detection on Chest X-Rays with Deep Learning](#):

Person or Device	F1	95% CI 2sigma	68% CI 1sigma
Radiologist 1	0.383	(0.309, 0.453)	(0.345, 0.417)
Radiologist 2	0.356	(0.282, 0.428)	(0.319, 0.392)
Radiologist 3	0.365	(0.291, 0.435)	(0.327, 0.399)
Radiologist 4	0.442	(0.390, 0.492)	(0.416, 0.467)
Radiologist Avg.	0.387	(0.330, 0.442)	(0.358, 0.414)
CheXNet	0.435	(0.387, 0.481)	(0.411, 0.458)
XR-PDA Max F1	0.408		

We do not calculate here the models 95% confidence interval for simplicity, and compare models statistical significance by assuming normal distribution and simply comparing the models F1 score to 1-sigma CIs calculated from 2-sigma ones. This model's score is higher and outside of 68% (1-sigma) CI for two radiologists out of four, which points to 'some statistical significance' of the given model.

Comparing the F1 scores themselves, this model achieves higher maximum F1 score than an average radiologist in the study. State of the art neural network, as well as one radiologist from the study, do achieve higher F1 score, but the model's performance is comparable and in many cases exceeds the performance of human radiologists (in terms of F1 score).

Furthermore, since the model does not have a high precision with any meaningful recall value, its usefulness tends to lie in its recall (and negative predictive value). Therefore, it makes sense to maximize recall and NPV even at the cost of small loss in precision. A good threshold value that achieves that is 0.377:

Device	F1	Precision	Sensitivity/Recall	Specificity	NPV
XR-PDA Max F1	0.408	0.290	0.689	0.579	0.856
XR-PDA T=0.377	0.404	0.268	0.818	0.441	0.905

If the model predicts negative, it is correct with 90.5% probability. If the model predicts positive, it is correct with 26.8% probability.

Out of all negative cases the model correctly classifies 44.1%, out of all positive cases it correctly classifies 81.8%.

In the presence of Infiltration disease, the algorithm works poorly because it is not able to differentiate between Pneumonia and infiltration.

In summary :

The algorithm performs with the following score :


- Max F1 score: 0.404. Outperform the radiologist and remember that F1 is test's accuracy the algorithm
- threshold: 0.377,
- precision: 0.268,
- recall: 0.818

The recall is very high which is great because few False Negative won't be detected. And precision is low so the opposite (more False Positive) But in the case of AI Healthcare domain it is better to say people are infected even if it is not than the opposite

4. Databases**Description of Training Dataset:**

Training dataset consisted of 2290 chest xray images, with a 50/50 split between positive and negative cases.

Example images:

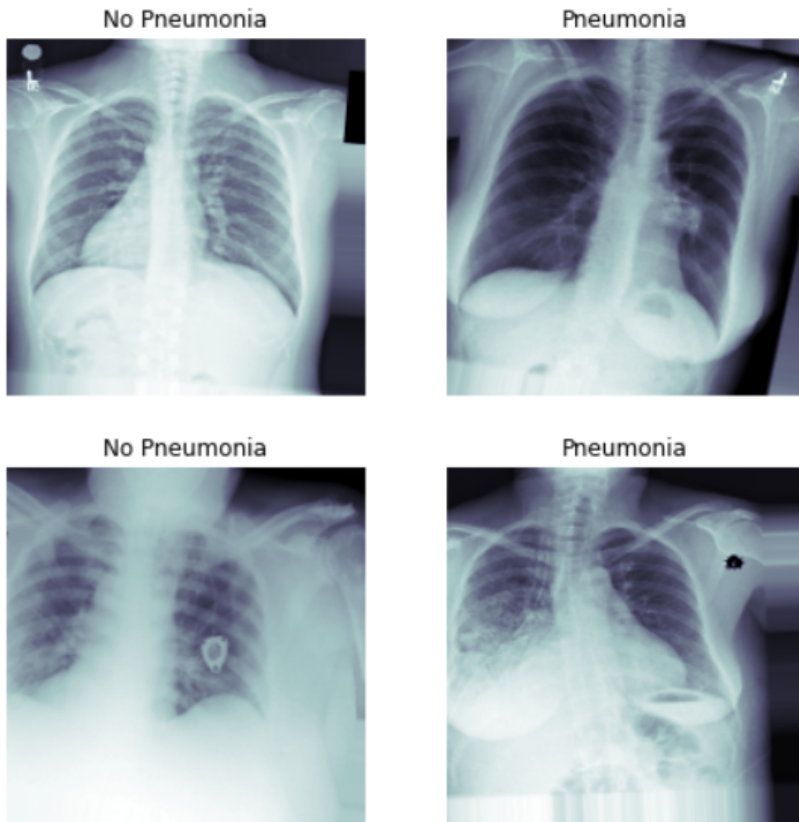

Description of Validation Dataset:

Validation dataset consisted of 1430 chest xray images, with 20/80 split between positive and negative cases, which more reflects the occurrence of pneumonia in the real world.

5. Ground Truth

The data is taken from a larger xray [dataset](#), with disease labels created using Natural Language Processing (NLP) mining the associated radiological reports. The labels include 14 common thoracic pathologies (Pneumonia being one of them):

- Atelectasis
- Consolidation
- Infiltration
- Pneumothorax
- Edema
- Emphysema
- Fibrosis
- Effusion
- Pneumonia
- Pleural thickening
- Cardiomegaly
- Nodule
- Mass
- Hernia

The biggest limitation of this dataset is that image labels were NLP-extracted so there could be some erroneous labels but the NLP labeling accuracy is estimated to be >90%.

The original radiology reports are not publicly available but more details on the labeling process can be found [here](#).

6. FDA Validation Plan

Patient Population Description for FDA Validation Dataset:

Ideal Dataset:

The ideal dataset would include Chest X-Rays of patient with Pneumonia but with an unbiased distribution with the following comorbid diseases (Consolidation, Edema, Effusion, Hernia). It will cover 50-50 Male/female and people mainly above 50 years old. No discrimination regarding the race.

And the following population subset is to be used for the FDA Validation Dataset:

- Both men and women
- Age 50 to 90
- Without known comorbid thoracic pathologies listed above

Ground Truth Acquisition Methodology:

The golden standard for obtaining ground truth would be to perform one of these tests (see this [Mayo Clinic Link](#)):

- Sputum test
- Pleural fluid culture

Yet, those tests are quite expensive, and in most cases diagnosis is concluded by the physician based on radiologist's analysis/description. Since the purpose of this device is assisting the radiologist (not replacing him), the ground truth for the FDA Validation Dataset can be obtained as an average of three practicing radiologists (as a widely used 'silver standard'). The same method is used in the mentioned paper.

Algorithm Performance Standard:

In terms of Clinical performance, the algorithm's performance can be measured by calculating F1 score against 'silver standard' ground truth as described above. The algorithm's F1 score should exceed **0.387** which is an average F1 score taken over three human radiologists, as given in [CheXNet: Radiologist-Level Pneumonia Detection on Chest X-Rayswith Deep Learning](#), where a similar method is used to compare device's F1 score to average F1 score over three radiologists.

A 95% confidence interval given in the paper for average F1 score of 0.387 is (0.330, 0.442), so algorithm's 2.5% and 97.5% percentiles should also be calculated to get 95% confidence interval. This interval, when subtracted the interval above for the average, should not contain 0, which will indicate statistical significance of its improvement of the average F1 score. The same method for assessing statistical significance is presented in the above paper.